

DABIT-012

Valium Reagent

www.fieldforensics.com



1) Remove Cap from Clear Tube and remove Sampler.


DO NOT TOUCH Sticky Sampler Tip.

2) **Dab Sampler Tip** on material to be tested. Dab firmly to collect a sample of suspect material.



3) Return Sampler to Clear Inner Tube, Sampler Tip first, and replace the cap.

4) Hold Sampler as shown with the **Tip Pointed Down**. **Squeeze Clear Tube AND Sampler** to break both ampoules. Squeeze to expel all of the liquid into the bottom of the tube. Agitate solution for 10 seconds to ensure sample tip comes into contact with solution.

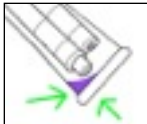

5) **Observe color of liquid.**

A **Purple** color indicates a **Positive** for **Valium**;
A light blue color is a negative test.

Observe for 60 seconds; if no **Purple** color then the test is negative.

Possible Irritant:

Avoid contact of liquids with eyes and mucous membranes.

November 2015 Rev. 8

FIELD FORENSICS INC.

DABIT-012

*Valium -
Valium/Ketamine/Rohypnol*

BZO
(Valium)

Control

DABIT-012: "Valium Reagent" - Tests for Valium, Ketamine and Rohypnol.

Full DABIT™ Color Chart Available at www.fieldforensics.com